


HeartSine® samaritan® PAD Trainer-Pak



Remove Electrode Pads

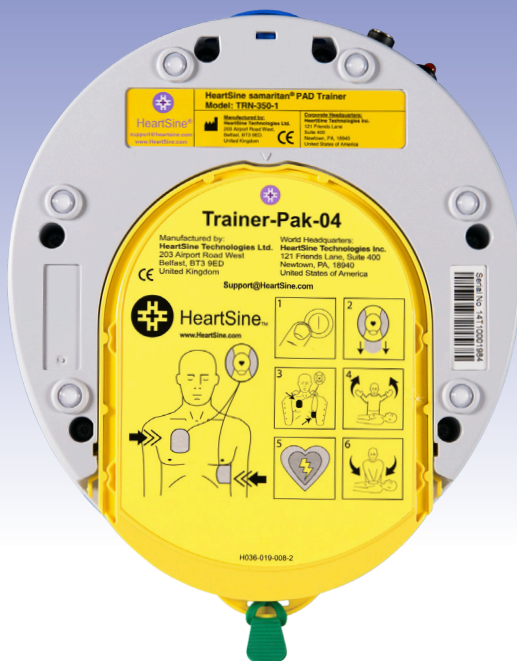


1. Pull green tab to open Trainer-Pak.




2. Pull pads from Trainer-Pak until the retractable cord lengthens.

Replace Electrode Pads

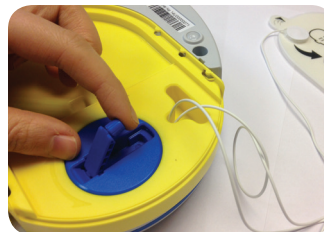

3. Lift blue rewind handle.

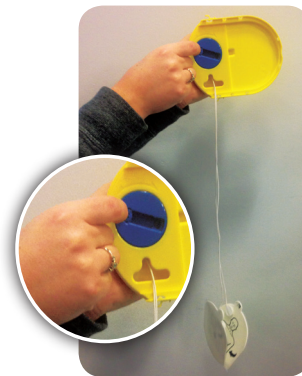

4. While holding the Trainer-Pak with the pad cord hanging, turn blue rewind handle clockwise to retract the cord.

Close Trainer-Pak



5. If the rewind mechanism locks up, eject the rewind mechanism to free the cord.



6. When cord is retracted, fold blue rewind handle into place.



7. Place pads in Trainer-Pak with peel tab at the top.



8. Slide cover over Trainer-Pak. Click cover into place.

EMEA/APAC

HeartSine Technologies, Ltd.
203 Airport Road West
Belfast, Northern Ireland
BT3 9ED
Tel: +44 28 9093 9400
Fax: +44 28 9093 9401
info@heartsine.com

**U.S./Americas**

HeartSine Technologies LLC
121 Friends Lane, Suite 400
Newtown, PA 18940
Toll Free: (866) 478 7463
Tel: +1 215 860 8100
Fax: +1 215 860 8192
info@heartsine.com

