

Accessories

HeartSine® samaritan® PAD Family

HeartSine samaritan PAD Accessories

Pad-Pak™

PAD-PAK-01

Adult Pad-Pak™ for patients > 8 years or 55 lbs/25 kg

PAD-PAK-02

Pediatric-Pak™ for patients 1 to 8 years or up to 55 lbs/25 kg

PAD-PAK-07

Adult Pad-Pak with TSO-C142a

HeartSine has taken an innovative approach to maintenance items by incorporating the AED batteries and electrodes into one single-use cartridge called the Pad-Pak. The advantage to the Pad-Pak is that there is only one date to monitor and only one item to replace on expiration or after use.

Each Pad-Pak has an industry-first four-year life, and provides 60 shocks or six hours of continuous monitoring. When battery power is down to less than 10 shocks, the samaritan PAD status indicator changes from green to red and emits an audible prompt to indicate the Pad-Pak should be replaced.

When used with the Pediatric-Pak, the samaritan PAD ensures the appropriate energy level is delivered for children ages 1-8 years old and weighing < 55 lbs/25 kg.

For use on airplanes, HeartSine offers the TSO-C142a certified Pad-Pak.



CPR Prep Kit

PAD-ACC-17

The CPR Prep Kit includes a CPR mask (15 mm port, one-way valve & filter), scissors, razor, two pair of non-latex exam gloves, antimicrobial wipe, towelette and plastic bag.



USB Data Cable

PAD-ACC-02

The custom USB Data Cable is for use with the Saver EVO™ data management software to download event data and to perform software upgrades to the samaritan PAD devices.



Saver EVO™ Data Management Software

Saver EVO is free icon-based, data management software for downloading, reviewing and printing ECG waveforms and event history data recorded in a samaritan PAD.



Features include:

- Comprehensive event documentation to facilitate physician oversight
- Utstein style reporting and report export capabilities
- User configuration setup and time synchronization
- Simple point and click operation

Saver EVO requires the optional HeartSine USB Data Cable to interface with the samaritan PAD. The Saver EVO software can be downloaded free of charge on the HeartSine website.

Carrying Case (Replacement)

PAD-BAG-01

Every HeartSine samaritan PAD is supplied with a custom-designed case with two zippered compartments. Each case has storage space to carry:

- HeartSine samaritan PAD AED
- CPR Prep Kit
- Additional Adult Pad-Pak or Pediatric-Pak
- Data Cable
- User Manual



Rigid Watertight Carrying Case

PAD-CAB-05

The new ultra-rugged carrying case, custom-fitted for the HeartSine samaritan PAD, features heavy-duty ABS construction that provides the ultimate protection and portability for your HeartSine AED. This long lasting dust and waterproof case (IP67) includes plenty of room for extra items, such as extra Pad-Paks and a first-aid kit.



Mobile AED Rescue Backpack

PAD-BAG-02

The HeartSine Mobile AED Rescue Backpack, custom made for the HeartSine samaritan Public Access Defibrillator (PAD), lets you create and transport a heart safe solution to virtually any location.

The HeartSine Mobile AED Rescue Backpack has room for everything you will need to keep members of your group heart safe. This customized backpack can readily carry a samaritan PAD AED with Pad-Pak, a Pediatric-Pak or extra Pad-Pak, CPR prep kit and first aid kit.




AED and first aid kit not included.

**Wall Cabinet with Alarm,
Custom Made for samaritan PAD**
PAD-CAB-04

Custom made for the samaritan PAD, this small profile wall cabinet sounds an alarm when opened.

Dimensions: 12"H x 11"W x 5"D

Weight: 5 lbs

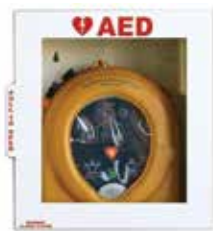

AED Wall Sign, Triangular
PAD-ACC-05

This three-dimensional wall sign is placed above the wall cabinet or wall bracket to call attention to the mounted samaritan PAD.



Wall Bracket

PAD-CAB-02

Made of metal, this premium wire wall bracket provides convenient storage for the samaritan PAD in a room or vehicle.



HeartSine samaritan PAD Trainer and Accessories

samaritan PAD Trainer

TRN-350-US (SAM 350P)

TRN-450-US (SAM 450P)

The samaritan PAD Trainer with remote control guides users through simulated analysis, energy delivery and prompted CPR with the look and feel of a live samaritan Public Access Defibrillator (PAD) without the actual charge and discharge of an electrical shock.

Includes Trainer-Pak with three pairs of training electrodes, battery charger and carrying case.



Remote (Replacement)

TRN-ACC-16 (SAM 350P)

TRN-ACC-19 (SAM 450P)

Multi-function remote control enables selection of among six AHA/ERC scenarios for live demonstrations or training. One remote control can even operate multiple trainer units at a time (within a range of 10-15 feet/3-5 meters).

The SAM 450P Trainer remote provides the training instructor with the flexibility to choose the appropriate CPR Rate Advisor prompts based on the student's compression rate.



Trainer-Pak (Replacement Electrode Cartridge)

TRAINER-PAK-04

With an easy-to-use rewind system, the Trainer-Pak Electrode Cartridge stores the electrode pads and cables. The handy rewind feature reduces time spent resetting the system.



Training Electrodes

10/pak TRN-ACC-02

25/pak TRN-ACC-03

Reusable, replaceable training electrode pads, for use on training manikins, are easy to place back on the pad liner.



Battery Charger (Replacement)

TRN-ACC-14

The Battery Charger charges the internal battery of the samaritan PAD Trainer in just four hours, which enables the Trainer to operate up to seven hours—a full day of training—on a single charge.



EMEA/APAC

HeartSine Technologies, Ltd.
203 Airport Road West
Belfast, Northern Ireland
BT3 9ED
Tel: +44 28 9093 9400
Fax: +44 28 9093 9401
info@heartsine.com



U.S./Americas

HeartSine Technologies LLC
121 Friends Lane, Suite 400
Newtown, PA 18940
Toll Free: (866) 478 7463
Tel: +1 215 860 8100
Fax: +1 215 860 8192
info@heartsine.com



The HeartSine products described in this brochure meet the European Medical Directive requirement.

UL Classified. See complete marking on product.

CAUTION: U.S. Federal law restricts this device to sale by or on the order of a licensed practitioner.

© 2016 HeartSine Technologies LLC. All rights reserved. H009-032-223-1

Lifesaving, Pure and Simple



www.heartsine.com

